



| | |
|------------------|--------------------------------------------|
| Product Name | Haming - The Future of Toothbrushes |
| Materials | Handle: Polypropylene Bristles: Nylon |
| Heat Resistance | 80 degrees Celsius |
| Contents | 1 toothbrush |
| Bristle Firmness | Regular |
| Price | 1,998 yen (tax included) |

Product Overview

This dentist-designed, revolutionary brush shape simultaneously polishes between the teeth and along the gumline.

Capable of simultaneously brushing between the teeth and along the gumline, "Haming" can also simultaneously brush dental brace brackets or wires from above and below. This revolutionary new dentist-designed brush shape provides superior brushing coverage. Effectively cleans along the gumline, an area of difficulty using traditional toothbrushes and interdental brushes. This revolutionary shape fits in between the teeth, and is also ideal for users of dental braces.

Features of dentist-designed "Haming"




Haming

A completely new idea in toothbrush shape


- One size suits any tooth, from infant teeth to permanent teeth.
- Easily brushes between teeth, as well as tooth surface and along the gumline, for thorough brushing coverage.
- Simultaneously brushes dental brace brackets and wires from above and below.
- For basic use, brush both top and bottom teeth holding the handle horizontally. With short side-to-side movements, add rotation and pressure as appropriate.

How to use

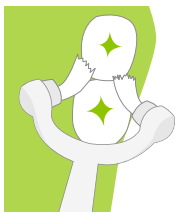


Rotating to the buccal (outward) side will simultaneously brush the occlusal (biting) surface and the buccal tooth surface.

Rotating to the buccal (outward) side will simultaneously brush the occlusal (biting) surface and the buccal tooth surface.

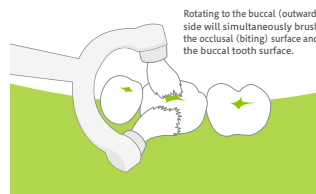


Rotating to the buccal (outward) side will simultaneously brush the occlusal (biting) surface and the buccal tooth surface.

Rotating to the buccal (outward) side will simultaneously brush the occlusal (biting) surface and the buccal tooth surface.

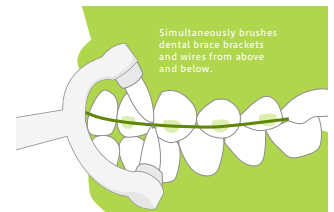

Simultaneously brushes dental brace brackets and wires from above and below.

Simultaneously brushes dental brace brackets and wires from above and below.

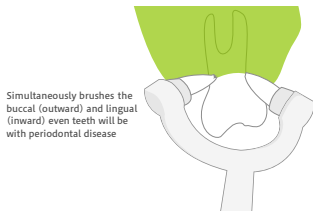

Simultaneously brushes the buccal (outward) and lingual (inward) even teeth will be with periodontal disease

Simultaneously brushes the buccal (outward) and lingual (inward) even teeth will be with periodontal disease



Brushing surface of the tongue with V-shaped brushes

Brushing surface of the tongue with V-shaped brushes

Contact Us about this Product

FREEDOM Co.,Ltd.

1st Floor, Yamazaki Building, 1-8-13 Nishi Kanagawa, Kanagawa Ward, Yokohama City, 221-0822

Tel: 045-548-3259 · Fax: 045-548-3269

E-mail: info@fdmjp.co.jp

WEB: <http://www.fdmjp.co.jp/en/>

Online Shop: <http://www.shopfdmjp.net/>

About
this Product

